

elephant HD Sensor



Patient Comfort


Rounded Edges: Designed for maximum patient comfort and

reduced anxiety.



Slim Profile: Enhances adaptability to oral structures, featuring a mere 4.5 mm thickness.



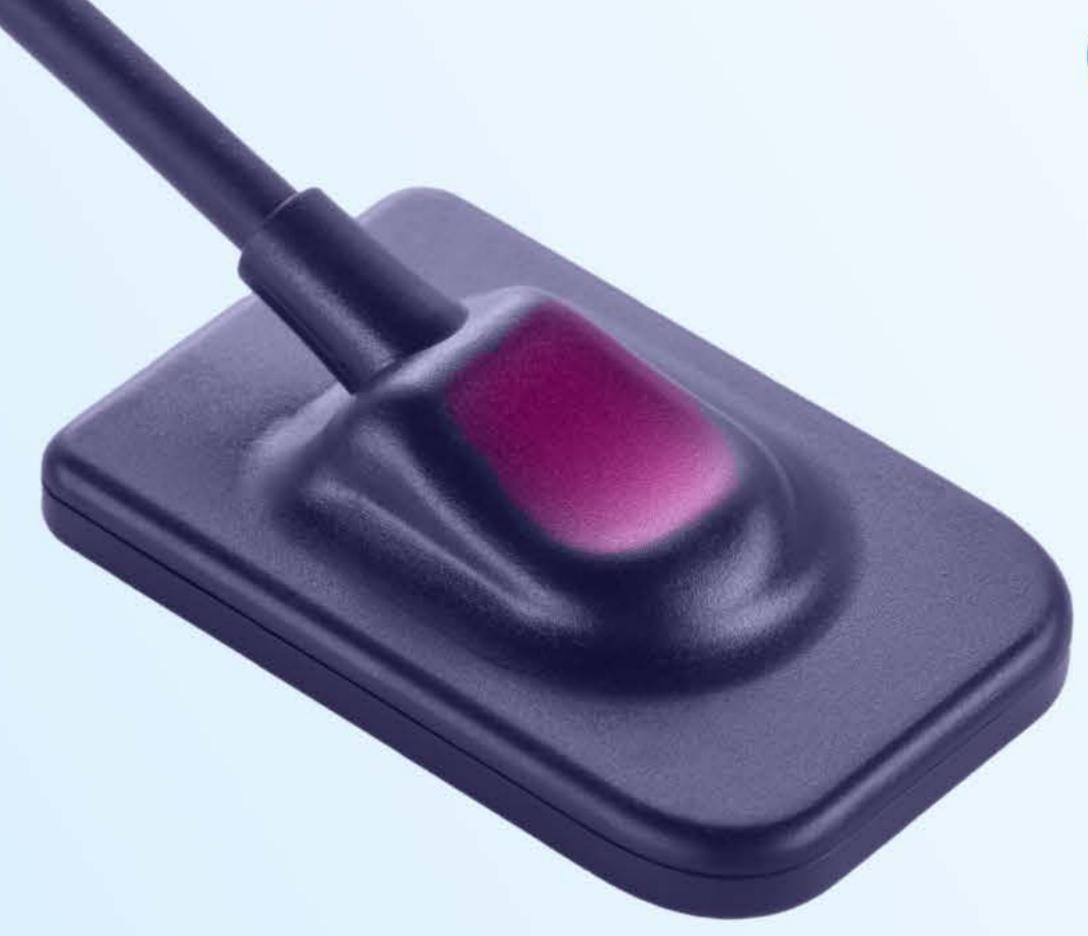

Concave Back Design

Ergonomically supports fingers for stable handling during procedures.

Connection Box

Positioned at the cable's end to minimize impact and protect against shocks.





Cable with High Reliability

Reinforced cord connection > 70,000 bending test (± 90°) passed

Sensor Sizes

Size 1.0

Not suitable in posterior region



Size 1.5
perfect balance of coverage and comfort



Size 2
Not comfortable due to its larger dimensions



Specifications	size 1	size 1.5	size 2
Dimensions	25.0*38.5*4.5	29.2*43.5*4.5	31.0*45.0*4.5
Active Area	20*30	24*33	26*36
Pixel	1000*1500	1200*1650	1300*1800
Sensor	APS CMOS	APS CMOS	APS CMOS
Grey Value	16Bit/0~65536	16Bit/0~65536	16Bit/0~65536
Power	5.0 V	5.0 V	5.0 V
Interface	USB 2.0	USB 2.0	USB 2.0

HD sensor specifications

Detector type: APS CMOS	Pixel size: 20um	
Scintillator. Csl	Grayscale: 16 bits	
Cable: replaceable	Resolution: 25lp/mm	
Line length: standard 3m	Protection level: 1p68	

- Hybrid ICs on circuit board
 Low noise integrated circuit design
- APS CMOS image sensor
 Low noise CMOS image sensors
 In-house design and process
- Csl scintillator

Directly deposited CsI Low Dose (LD) and High Resolution (HR) options

Robust and safety design

Biocompatible IP68 protection Easy infection control

Encapsulation layer

Protects scintillator



The highest level of protection against water and dust in the world.

Degree of Protection against SOLID object.

O No protection




solid object >50mm as a hand




solid object>12.5mm as a finger



solid object >2.5mm as a wire

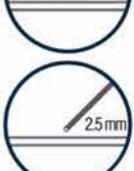


solid object >1mm as a strap

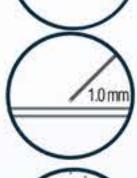

5 dust protection



6 no ingress of dust







Dust sight.

No ingress of

Protected against the affects of permanent submersion in water (up to 13 feet)

Degree of Protection against LIQUID (Water).

O No protection



water drops



2 water drops at 15°



3 water spray at 60°




water splashing




5 water jets



6 powerful jet and heavy seas

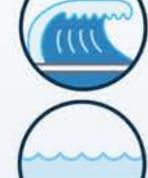

7 Temporary submersion (30min. at 3 feet)

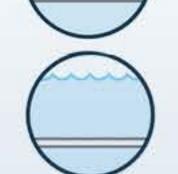

9 permanent submersion (up to 13 feet)



Brillant
Waterproof
performance



Sensor Software

User-Friendly Interface

Designed to be intuitive and easy to navigate, ensuring a smooth user experience for dentists.


Simplified Operations


Provides a straightforward approach to image capture and analysis, saving time and effort.



Enhanced Performance

Offers features that optimize diagnostics and improve treatment planning, supporting accurate and effective dental care.

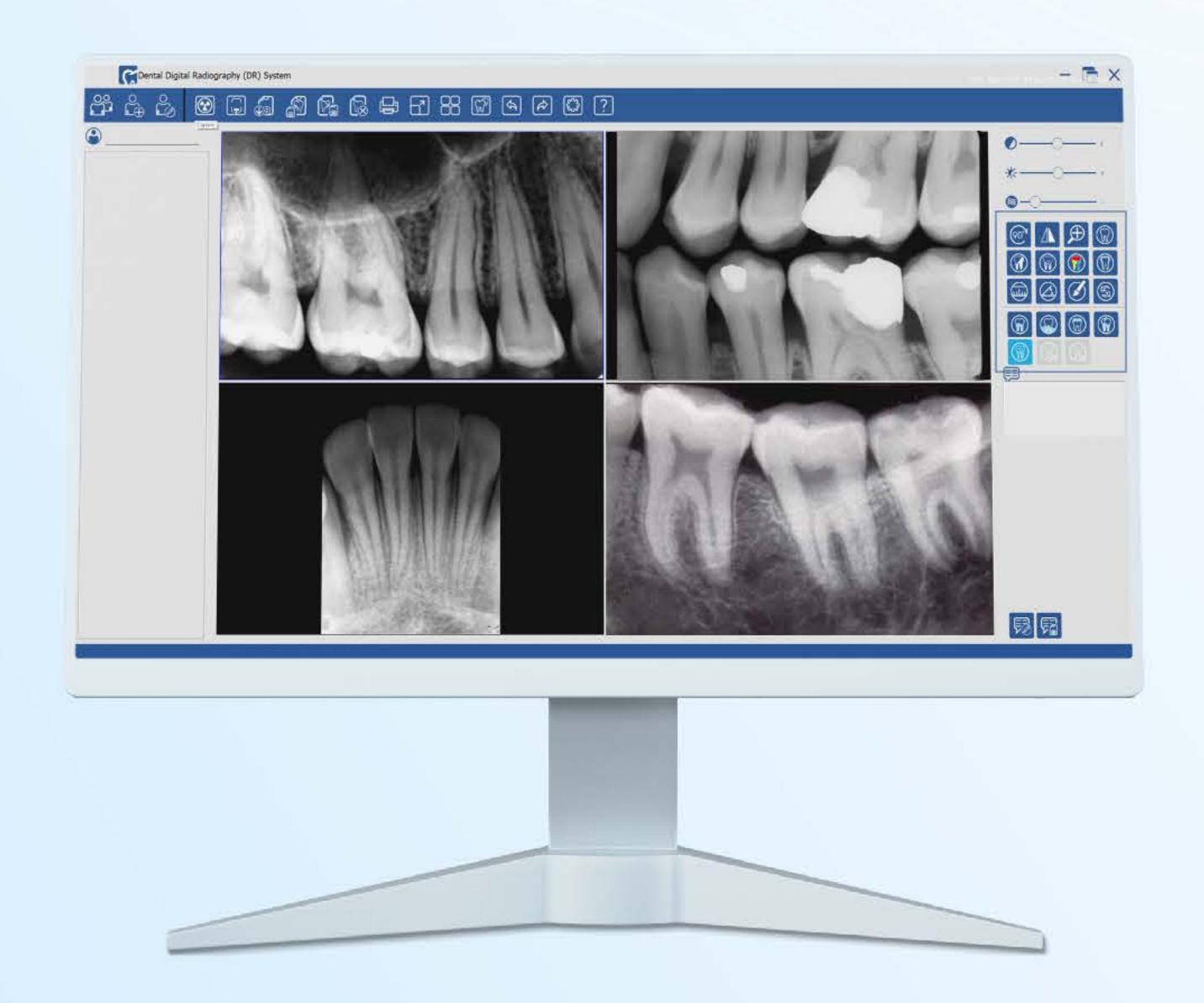

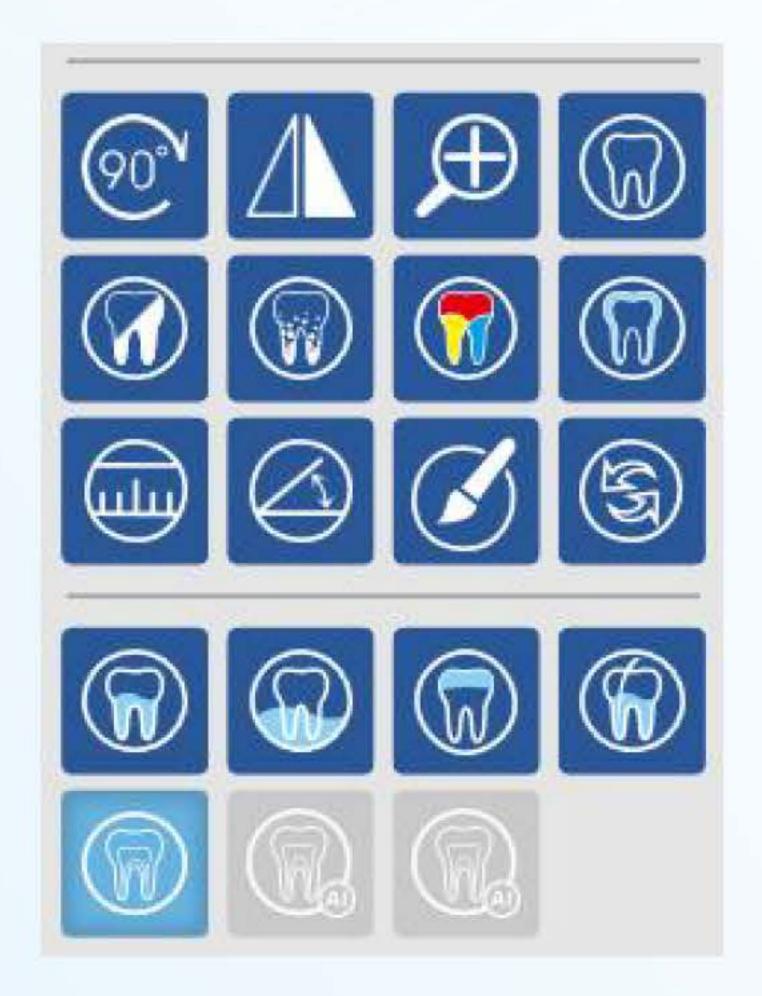

Icons for the ability to edit and customize the settings of the captured images to suit your needs.



CONTACT US

- 0085292083786 00852440568345
- info@elephant-dent.com
- Address:
 ROOM 1805 BLK 118 F EVERGAIN PLAZA
 88 CONTAINER PORT ROAD, HONG KONG

